

Chiral Separation of Nadolol on Lux 5µm i-Amylose-1 in NP


Column: Lux® 5 µm i-Amylose-1, LC Column 250 x 4.6 mm, Ea
Dimensions: 250 x 4.6 mm ID
Order No: 00G-4762-E0
Elution Type: Isocratic
Eluent A: 0.1% DEA in Hexane
Eluent B: 0.1% DEA in IPA
Gradient

Step No.	Time (min)	Pct A	Pct B
1	0	85	15

Profile:
Flow Rate: 1 mL/min
Col. Temp.: ambient
Detection: Electrospray Mass Spec (ESMS) @ 0.0000000000 (ambient)



Products used in this application:

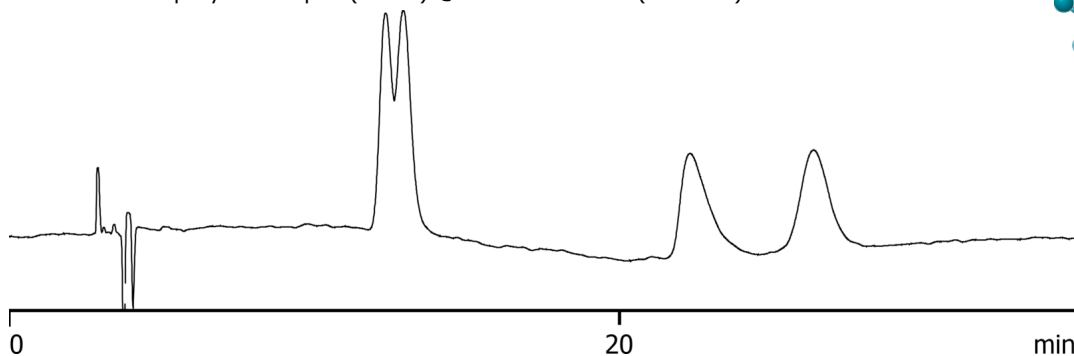

ANALYTES:

- 1 Nadolol
Retention Time: 12.355 min
- 2 Nadolol
Retention Time: 12.915 min
- 3 Nadolol
Retention Time: 22.312 min
- 4 Nadolol
Retention Time: 26.359 min

