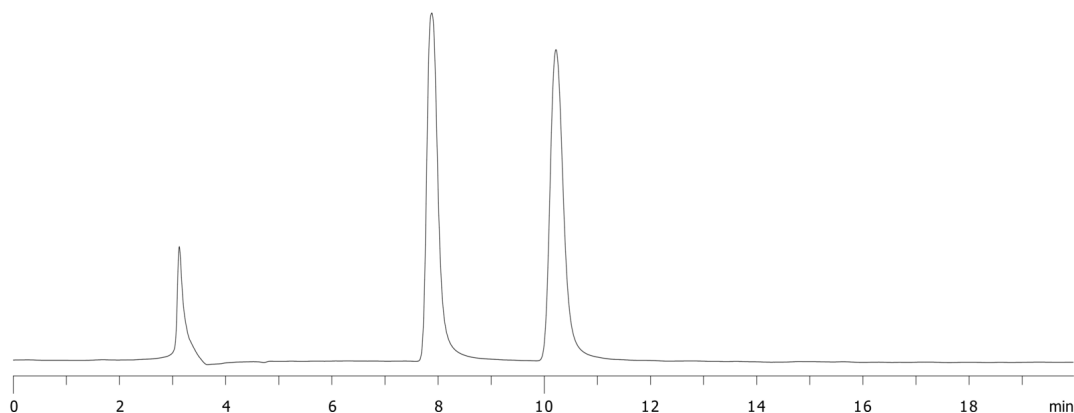


Chiral Separation of Lorazepam on Lux 5 μ m i-Amylose-1 250 x 4.6mm in NP**Column:** Lux[®] 5 μ m i-Amylose-1, LC Column 250 x 4.6 mm, Ea**Dimensions:** 250 x 4.6 mm ID**Order No:** 00G-4762-E0**Elution Type:** Isocratic**Eluent A:** n-Hexane**Eluent B:** Acetone**Gradient**

Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	60	40

Flow Rate: 1 mL/min**Col. Temp.:** 25 °C**Detection:** Electrospray Mass Spec (ESMS) @ 0.000000000 (ambient)

24453

**Products used in this application:****ANALYTES:****1** Lorazepam

Retention Time: 7.6 min

2 Lorazepam

Retention Time: 9.9 min

