

NSAIDs extracted on Strata DE 200uL 96 Well

Column: Kinetex® 2.6 µm C18 100 Å, LC Column 50 x 2.1 mm, Ea

Dimensions: 50 x 2.1 mm ID

Order No: 00B-4462-AN

Elution Type: Gradient

Eluent A: 0.1% formic acid in water

Eluent B: 0.1% formic acid in acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	90	10
	2	2.5	0	100
	3	3.5	0	100
	4	3.51	90	10
	5	5.51	90	10


Flow Rate: 500

Col. Temp.: ambient

Detection: Electrospray Mass Spec (ESMS) @ 0.0000000000 (ambient)
App ID 2481



Products used in this application:

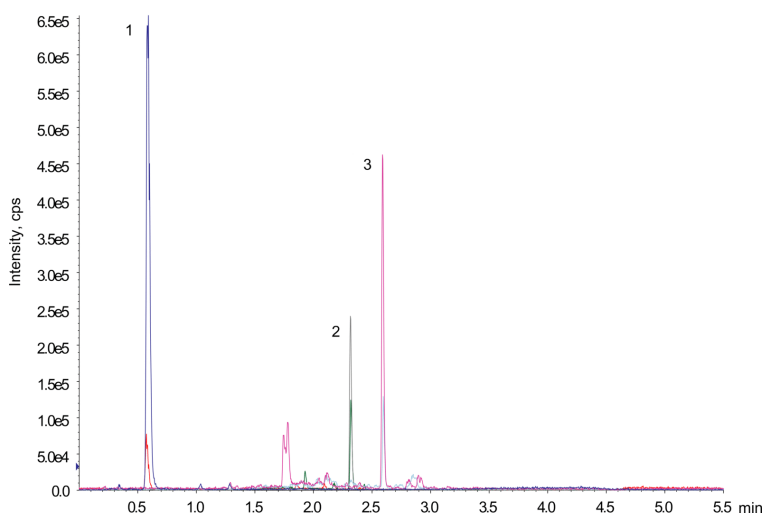

ANALYTES:

- 1** acetaminophen
Retention Time: 0.59 min
- 2** acetaminophen D7
Retention Time: 0.59 min
- 3** naproxen
Retention Time: 2.32 min
- 4** naproxen D3
Retention Time: 2.32 min
- 5** diclofenac
Retention Time: 2.6 min
- 6** diclofenac D4
Retention Time: 2.6 min

