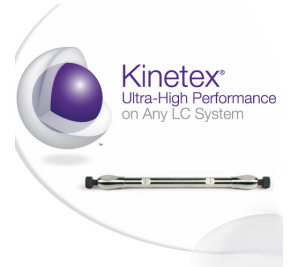
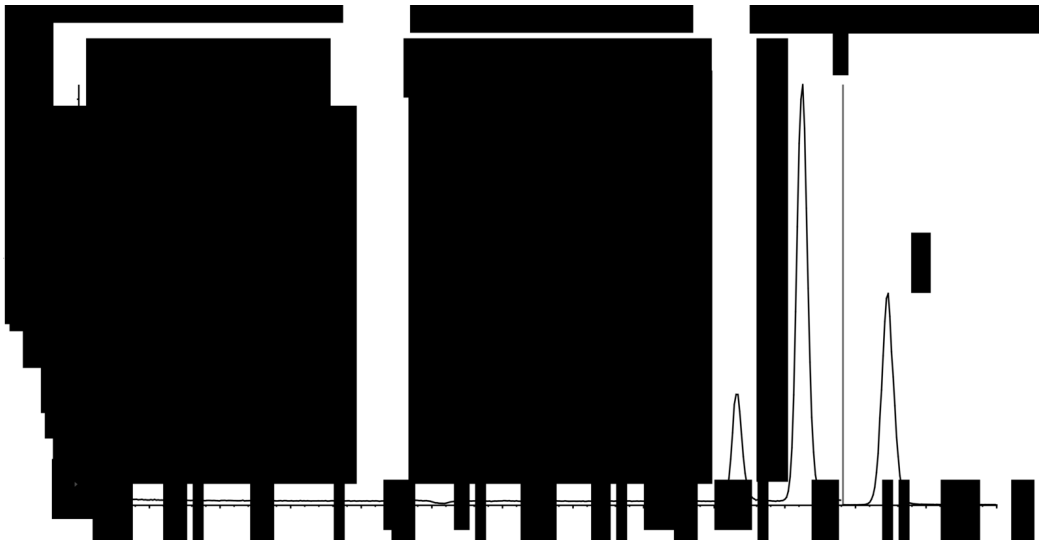


Analysis of Uracil and metabolites on Kinetex 2.6 µm PS C18 150x3.0 Column**Column:** Kinetex® 2.6 µm PS C18 100 Å, LC Column 150 x 3 mm, Ea**Dimensions:** 150 x 3 mm ID**Order No:** 00F-4780-Y0**Elution Type:** Gradient**Eluent A:** 0.1% Formic acid in Water**Eluent B:** Methanol

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	93	7
	2	6	93	7
	3	6.01	50	50
	4	7	50	50
	5	7.01	93	7
	6	12	93	7

Flow Rate: 0.2 mL/min**Col. Temp.:** 15 °C**Detection:** LC/MS/MS @ 0.0000000000 m/z (mass per charge) (15 °C)**Analyst Note:** Analyte diluted with initial mobile phase**Products used in this application:****ANALYTES:**

- 1 5,6-dihydrouracil (UH2)

Retention Time: 4.6 min

- 2 Uracil

Retention Time: 5.1 min

- 3 5-fluorouracil

Retention Time: 5.7 min

