

GC-MS/MS Separation of Hexachlorodibenzodioxins (HxCDD) on ZB-Dioxin GC column with 5m Zebron Guardian

Column: Zebron™ ZB-Dioxin w/GUARDIAN™, 5 m. GC Cap. Column 60 m x 0.25 mm x 0.20 µm, E

Phase:

Dimensions: 60 meters x 0.25 mm x 0.2 µm


Order No: 7KG-G045-10-GGA

Oven Profile: to 220°C @ 30°C/min 220°C for 15.5 min to 240°C @ 2°C/min 240°C for 4.5 min to 270°C @ 5°C/min 270°C for 5 min to 330°C @ 15°C/min


Carrier Gas: Constant Flow Helium, 1.056 mL/min

Injection: Splitless :1 1 µL @ 330°C

Detection: GC/MS/MS (340°C)

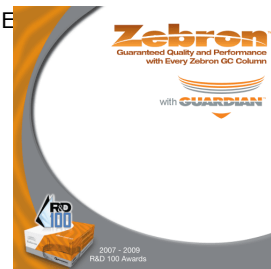

Products used in this application:

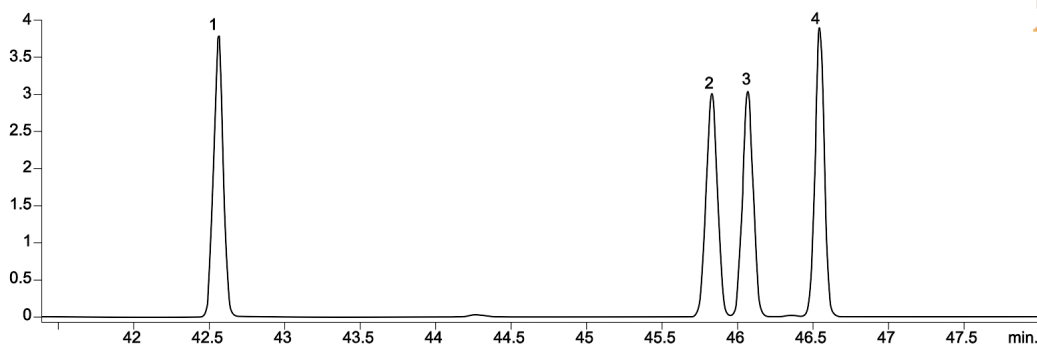**ANALYTES:**

1 1,2,4,6,7,9-HxCDD

Retention Time: 42.5 min

2 1,2,3,4,7,8-HxCDD

Retention Time: 45.8 min

3 1,2,6,7,8-HxCDD

Retention Time: 46.2 min

4 1,2,3,7,8,9-HxCDD

Retention Time: 46.5 min

