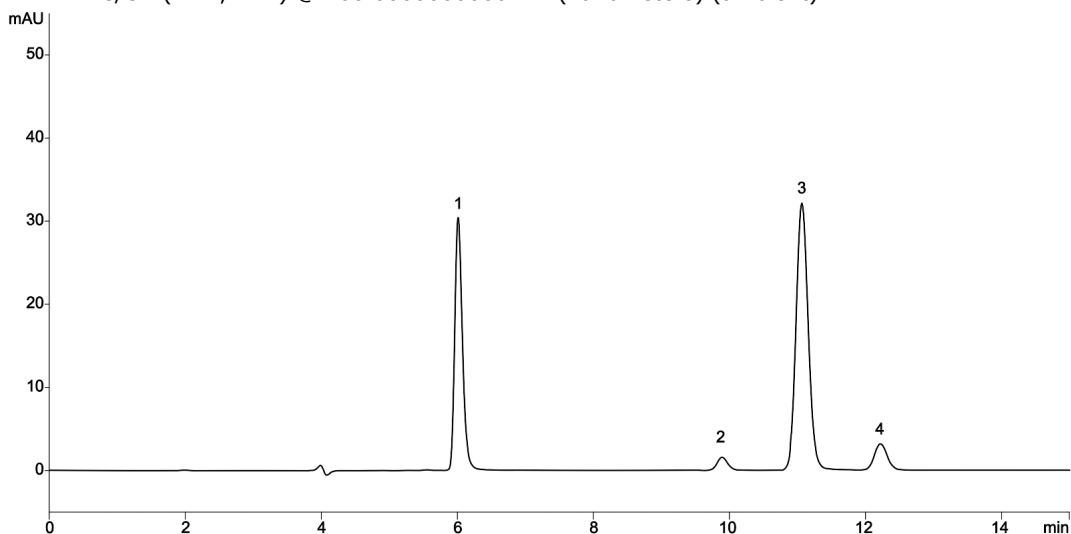


Fluconazole Standard solution (mix of 4) USP on Inertsil 3 μ m ODS-3 150x4.6mm**Column:** N/A**Dimensions:** N/A**Order No:** N/A**Elution Type:** N/A**Eluent A:** Water:Acetonitrile 80:20 (Pre-mixed)**Gradient** **Step No.** **Time (min)** **Pct A****Profile:** **1** 15 100**Flow Rate:** 0.5 mL/min**Col. Temp.:** 40 °C**Detection:** LC/UV (DAD, PDA) @ 260.000000000 nm (nanometers) (ambient)**ANALYTES:****1** Fluconazole Related Compound A

Retention Time: 6.013 min

2 Fluconazole Related Compound B

Retention Time: 9.901 min

3 Fluconazole Related Compound C

Retention Time: 11.065 min

4 Fluconazole

Retention Time: 12.222 min

