

HPLC Application

ID No.: 27718

PhEur Atazanavir Sulfate - Impurity K; Reference solution D on Luna Omega 3u C18 150x4.6mm using Agilent 12

Column: Luna® Omega 3 µm C18 100 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4784-E0

Elution Type: Gradient

Eluent A: mix 25 volumes of acetonitrile and 75 volumes of a freshly prepared 2.73 g/L solution of


Eluent B: mix 25 volumes of a freshly prepared 2.73 g/L solution of potassium dihydrogen phosphate p

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	5.6	95	5
	3	8.6	0	100
	4	14.6	0	100
	5	16	95	5
	6	23	95	5


Flow Rate: 1 mL/min

Col. Temp.: 25 °C

Detection: LC/UV (DAD, PDA) @ 215.000000000 nm (nanometers) (ambient)

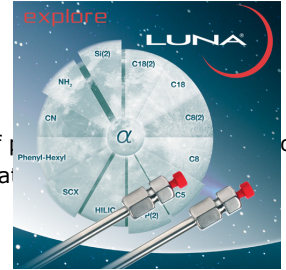

Products used in this application:

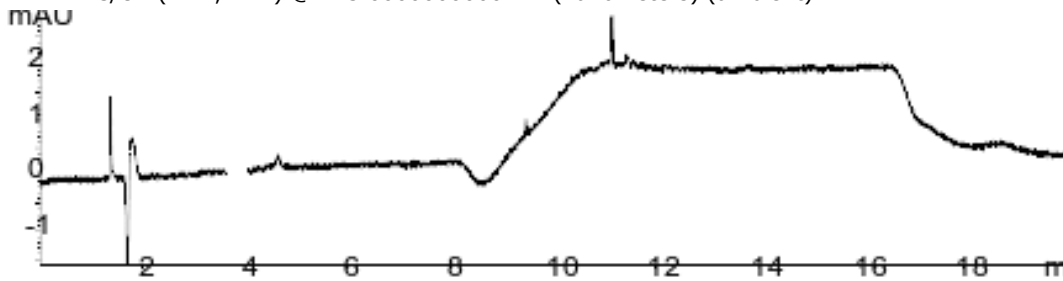

ANALYTES:

1 Impurity K

Retention Time: 3.87 min

