

Prednisolone - USP Acceptable Modified Method

Column: Luna® 5 µm C8(2) 100 Å, LC Column 100 x 4.6 mm, Ea

Dimensions: 100 x 4.6 mm ID

Order No: 00D-4249-E0

Elution Type: Isocratic

Eluent A: Butyl chloride:water-saturated butyl chloride:tetrahydrofuran:methanol:glacial acetic acid

Gradient	Step No.	Time (min)	Pct A
Profile:	1	0	100

Flow Rate: 1.5 mL/min

Col. Temp.: ambient

Detection: Electrospray Mass Spec (ESMS) @ 0.0000000000 (ambient)


Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C8 4 x 3.0mm, 10/Pk Part No.: AJ0-4290

- Holder Part No.: KJ0-4282



Products used in this application:

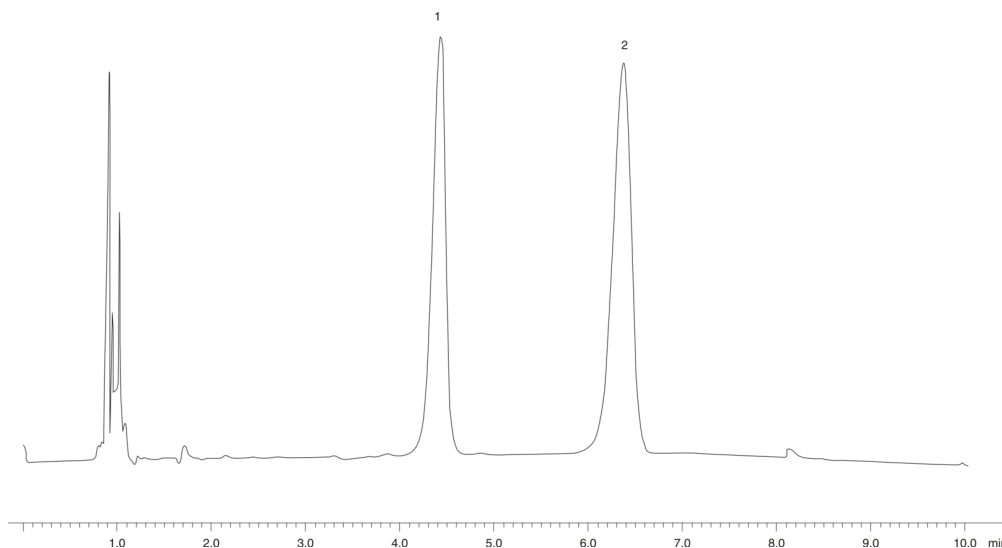

ANALYTES:

1 Betamethasone

2 Prednisolone

