

## Amitriptyline HCl - Alternative Method

**Column:** Luna® 5 µm Phenyl-Hexyl 100 Å, LC Column 150 x 4.6 mm, Ea

**Dimensions:** 150 x 4.6 mm ID

**Order No:** 00F-4257-E0

**Elution Type:** Isocratic

**Eluent A:** 20mM Ammonium acetate pH 4.9

**Eluent B:** Acetonitrile

**Eluent C:** Methanol

| Gradient Profile: | Step No. | Time (min) | Pct A | Pct B | Pct C |
|-------------------|----------|------------|-------|-------|-------|
| 1                 | 0        | 40         | 10    | 50    |       |

**Flow Rate:** 2 mL/min

**Col. Temp.:** ambient

**Detection:** Electrospray Mass Spec (ESMS) @ 0.0000000000 (ambient)


**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Phenyl 4 x 3.0mm, 10/Pk Part No.: AJ0-4351

- Holder Part No.: KJ0-4282



Products used in this application:

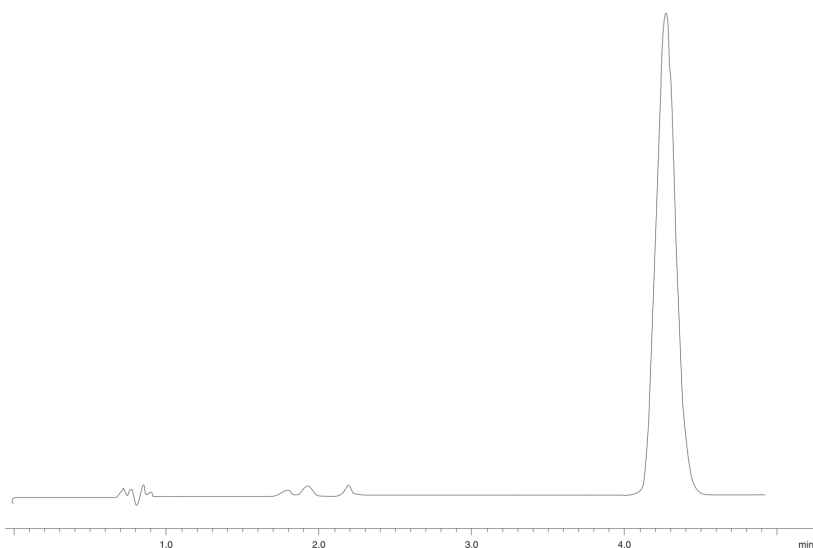

### ANALYTES:

1 Amitriptyline

