


**EPA Method 624 Calibration Mix****Column:** Zebron<sup>™</sup> ZB-624, GC Cap. Column 60 m x 0.25 mm x 1.40  $\mu$ m, Ea**Phase:** 6% Cyanopropylphenyl 94% Dimethylpolysiloxane**Dimensions:** 60 meters x 0.25 mm x 1.4  $\mu$ m**Order No:** 7KG-G005-27**Oven Profile:** 40 $^{\circ}$ C for 4 min 160 $^{\circ}$ C for 28 min**Carrier Gas:** Constant Flow Hydrogen, 1.03 mL/min**Injection:** Split 100.0000000000:1 0  $\mu$ L @ 250 $^{\circ}$ C**Detection:** Refractive Index (255 $^{\circ}$ C)**Products used in this application:**

## EPA Method 624 Calibration Mix

### ANALYTES:

- 1 1,1-Dichloroethene
- 2 Methylene chloride
- 3 1,1-Dichloroethane
- 4 Carbon tetrachloride
- 5 Benzene
- 6 Trichloroethene
- 7 1,2-Dichloropropane
- 8 2-Chlorethyl vinyl ether
- 9 1,1,2-Trichloroethane
- 10 Tetrachloroethene
- 11 Dibromochloromethane
- 12 Chlorobenzene

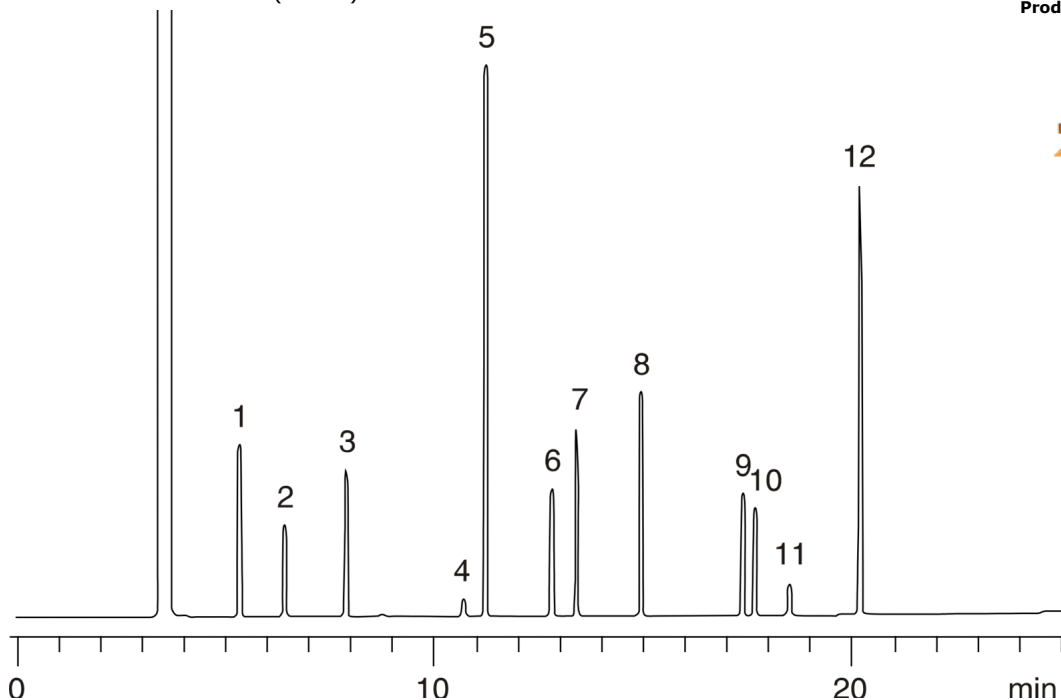