


EPA Method 608/CLP Pesticides**Column:** Zebron™ ZB-1701, GC Cap. Column 15 m x 0.25 mm x 0.25 µm, Ea**Phase:** 14% Cyanopropylphenyl 86% Dimethylpolysiloxane**Dimensions:** 15 meters x 0.25 mm x 0.25 µm**Order No:** 7EG-G006-11**Oven Profile:** 200Å°C for 5 min 250Å°C for 15 min**Carrier Gas:** Constant Flow Hydrogen, 1.47 mL/min**Injection:** Split 100.0000000000:1 0 µL @ 250°C**Detection:** Refractive Index (285°C)**Products used in this application:**

EPA Method 608/CLP Pesticides

ANALYTES:

- 1** alpha-BHC
- 2** gamma-BHC
- 3** Heptachlor
- 4** Aldrin
- 5** beta-BHC
- 6** delta-BHC
- 7** Heptachlor epoxide
- 8** Endosulfan I
- 9** 4,4'-DDE
- 10** Dieldrin
- 11** Endrin
- 12** 4,4-DDD
- 13** Endosulfan II
- 14** 4,4'-DDT
- 15** Endrin aldehyde
- 16** Endosulfan sulfate

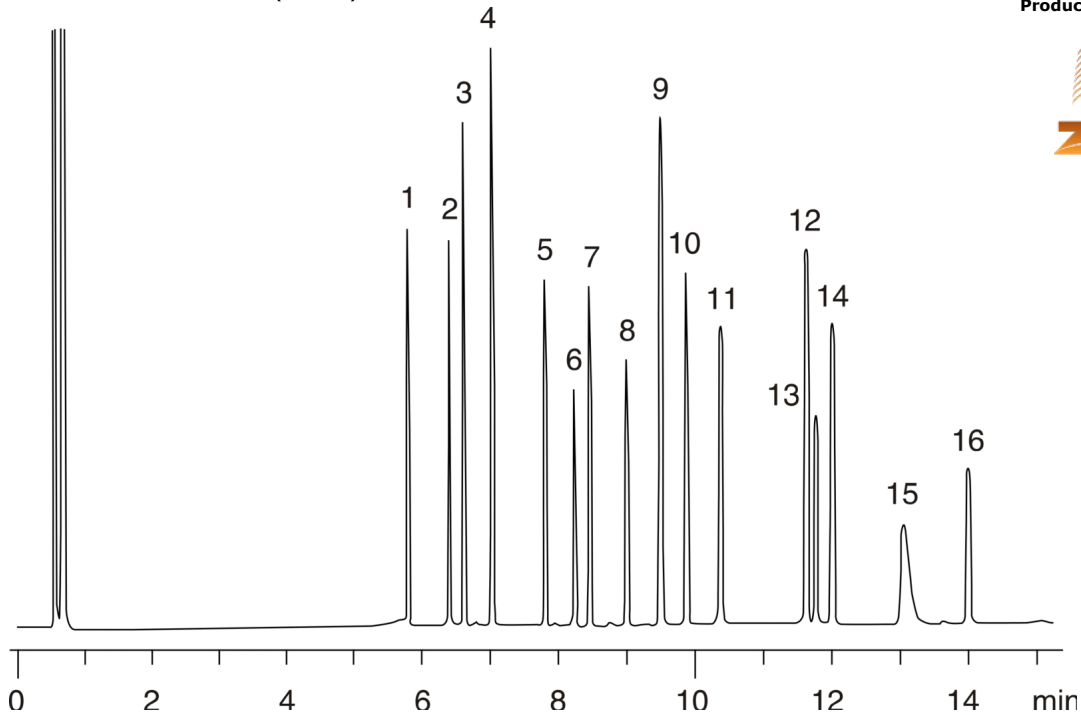