

Sugars in Honey on Rezex ROA

Column: Rezex™ ROA-Organic Acid H+ (8%), LC Column 300 x 7.8 mm, Ea

Dimensions: 300 x 7.8 mm ID

Order No: 00H-0138-K0

Elution Type: Isocratic

Eluent A: 0.01 N H₂SO₄

Gradient Profile:	Step No.	Time (min)	Pct A
1	1	0	100

Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: Differential refractometer @ 0.0000000000 (ambient)

Analyst Note: Standard solution (A) vs. Honey sample (B)

SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Carbo-H 4 x 3.0mm, 10/Pk Part No.: AJ0-4490


- Holder Part No.: KJ0-4282

REZEX™
HPLC COLUMNS for Carbohydrates, Oligosaccharides, and Organic Acid Analysis



Products used in this application:

REZEX™

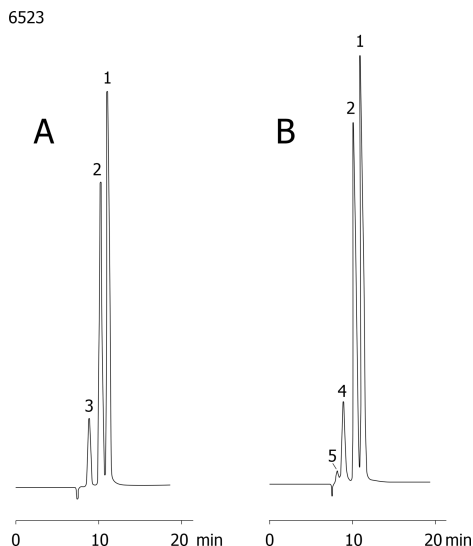

ANALYTES:

- 1 Fructose
- 2 Glucose
- 3 Sucrose
- 4 Disaccharide
- 5 Trisaccharide

