


Free Fatty Acids**Column:** Zebtron™ ZB-FFAP, GC Cap. Column 30 m x 0.32 mm x 0.25 µm, Ea**Phase:** Nitroterephthalic Acid Modified Polyethylene Glycol**Dimensions:** 30 meters x 0.32 mm x 0.25 µm**Order No:** 7HM-G009-11**Oven Profile:** 100Å°C for 2 min 180Å°C for 12 min 180Å°C for 16 min**Carrier Gas:** Constant Flow Helium, 1.69 mL/min**Injection:** Split 100.000000000:1 0 µL @ 250°C**Detection:** Refractive Index (350°C)

8186

**Products used in this application:****ANALYTES:**

- 1 Acetic acid
- 2 Propionic acid
- 3 Isobutyric acid
- 4 n-Butyric acid
- 5 Isovaleric acid
- 6 n-Valeric acid
- 7 Isocaproic acid
- 8 n-Caproic acid
- 9 n-Heptanoic acid

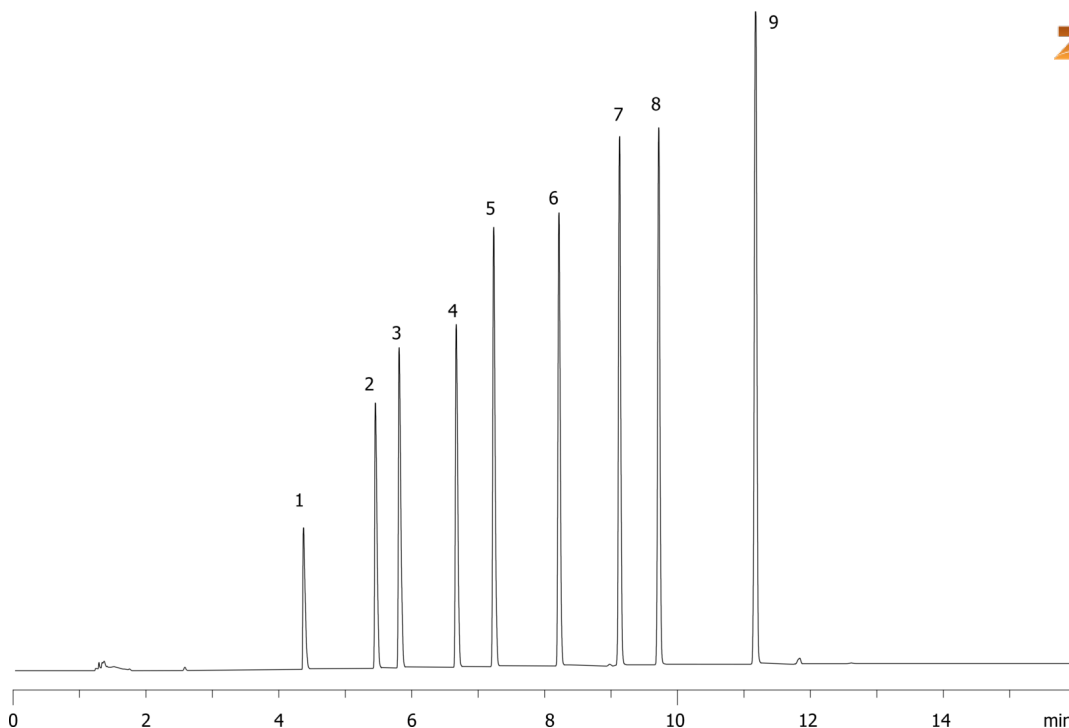