


Organic Acids**Column:** Zebtron™ ZB-FFAP, GC Cap. Column 30 m x 0.32 mm x 0.25 µm, Ea**Phase:** Nitroterephthalic Acid Modified Polyethylene Glycol**Dimensions:** 30 meters x 0.32 mm x 0.25 µm**Order No:** 7HM-G009-11**Oven Profile:** 100Å°C for 2 min 250Å°C for 20.75 min 250Å°C for 50.75 min**Carrier Gas:** Constant Flow Helium, 1.69 mL/min**Injection:** Split 20.0000000000:1 0 µL @ 250°C**Detection:** Refractive Index (350°C)**Products used in this application:**

Organic Acids

ANALYTES:

- 1 Acetic acid
- 2 Propionic acid
- 3 Isobutyric acid
- 4 n-Butyric acid
- 5 Isovaleric acid
- 6 n-Valeric acid
- 7 Isocaproic acid
- 8 n-Caproic acid
- 9 n-Heptanoic acid
- 10 n-Nonanoic acid
- 11 n-Undecanoic acid
- 12 n-Tridecanoic acid
- 13 n-Pentadecanoic acid (C15)
- 14 n-Heptadecanoic acid (C17)
- 15 n-Nonadecanoic acid (C19)
- 16 n-Heneicosanoic acid (C21)
- 17 n-Tricosanoic acid (C23)

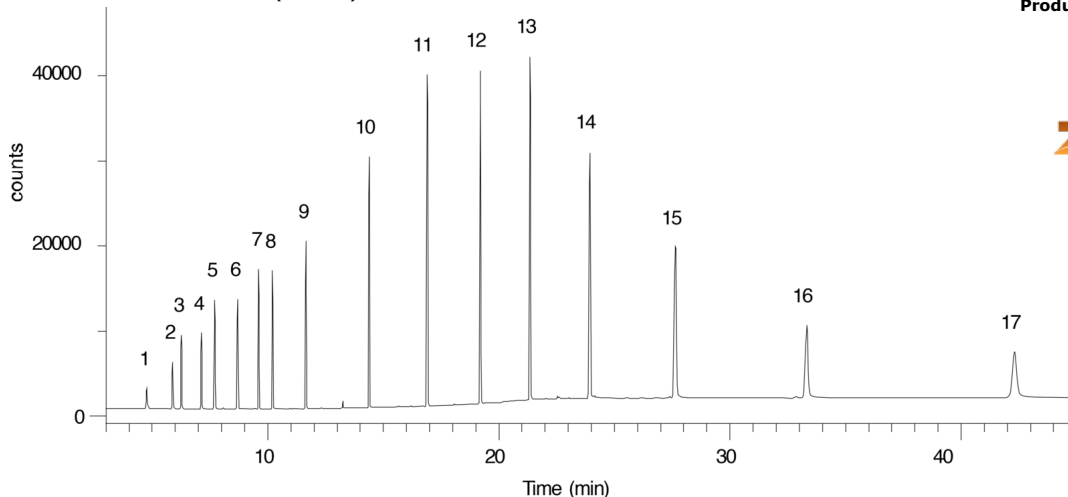