

Dopamine HCl Injection - USP Acceptable Modified Method

Column: Luna® 5 µm C18(2) 100 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4252-E0

Elution Type: Isocratic

Eluent A: Water with 1% Acetic acid and 5mM 1-octane sulfonate

Eluent B: Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
1	0	87	13	

Flow Rate: 2 mL/min

Col. Temp.: ambient

Detection: Differential refractometer @ 0.0000000000 (ambient)


Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C18 4 x 3.0mm, 10/Pk Part No.: AJ0-4287

- Holder Part No.: KJ0-4282



Products used in this application:

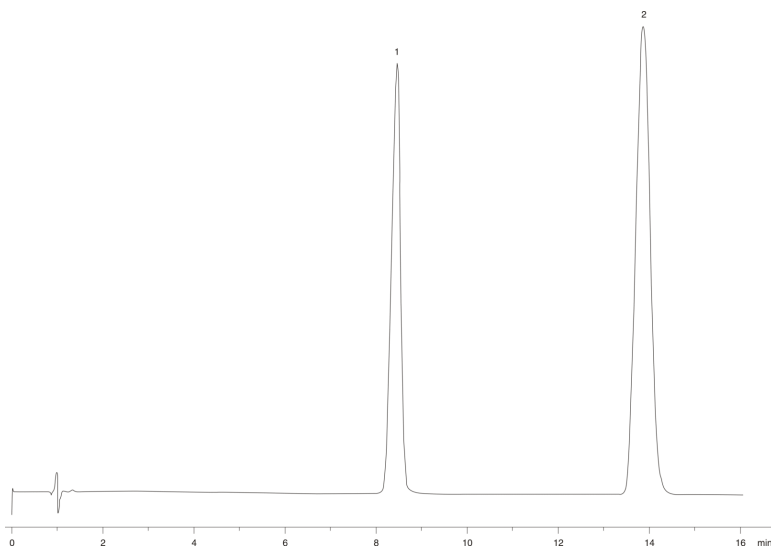

ANALYTES:

- 1 Benzoic acid
- 2 Dopamine

