

## Chemotactic Peptides

**Column:** Jupiter® 5 µm C5 300 Å, LC Column 250 x 4.6 mm, Ea

**Dimensions:** 250 x 4.6 mm ID

**Order No:** 00G-4052-E0

**Elution Type:** Gradient

**Eluent A:** 0.1% TFA in water

**Eluent B:** 0.1% TFA in acetonitrile

| Gradient Profile: | Step No. | Time (min) | Pct A | Pct B |
|-------------------|----------|------------|-------|-------|
|                   | 1        | 0          | 60    | 40    |
|                   | 2        | 20         | 20    | 80    |

**Flow Rate:** 1 mL/min

**Col. Temp.:** ambient

**Detection:** Differential refractometer @ 0.0000000000 (ambient)


**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Widespore C5 4 x 3.0mm, 10/Pk Part No.: AJ0-4327

- Holder Part No.: KJ0-4282



Products used in this application:

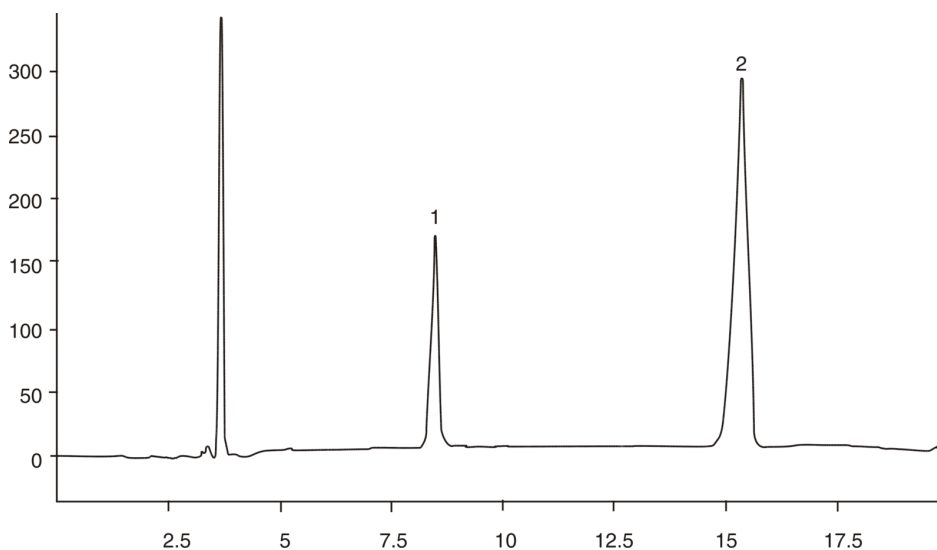

### ANALYTES:

1 N-Formyl-Nle-Leu-Phe-Nle-Tyr-Lys

2 Icararia

