

Hydrocortisone Acetate - USP Method

Column: Luna® 10 µm Silica (2) 100 Å, LC Column 250 x 4.6 mm, Ea
Dimensions: 250 x 4.6 mm ID
Order No: 00G-4091-E0
Elution Type: Isocratic
Eluent A: Butyl chloride:water-saturated butyl chloride:tetrahydrofuran:methanol:glacial acetic acid
Gradient


Step No.	Time (min)	Pct A
1	0	100

Profile:
Flow Rate: 1 mL/min
Col. Temp.: ambient
Detection: Differential refractometer @ 0.0000000000 (ambient)
Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Silica 4 x 3.0mm, 10/Pk Part No.: AJ0-4348
- Holder Part No.: KJ0-4282



Products used in this application:

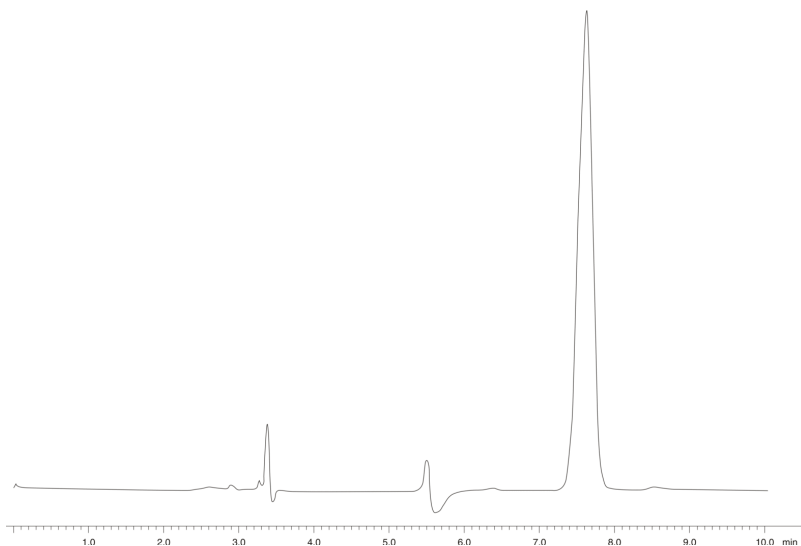

ANALYTES:

- Hydrocortisone acetate
Retention Time: 7.589 min

